

Q-Plex™ Mouse Cytokine Screen HS Stripwells (16-Plex)

The Q-Plex™ Mouse Cytokine HS Stripwells (16-Plex) is a fully quantitative ELISA-based chemiluminescent assay allowing the concurrent measurement of sixteen biomarkers or analytes: GMCSF, IFN γ , IL-1 α , IL-1 β , IL-2, IL-3, IL-4, IL-5, IL-6, IL-10, IL-12p70, IL-17a, MCP-1, MIP-1 α , RANTES, and TNF α .

Each kit contains a 96-well plate, featuring the relevant biomarker panel, and all the reagents required to perform testing. Q-Plex™ plates are built by absorbing sixteen distinct capture antibodies in a defined array to the bottom of each well. Our high quality reagents help ensure the accuracy of your results.

Catalog # Product

112149MS	Q-Plex™ Mouse Cytokine Screen HS Stripwells (16-Plex)	GMCSF	IFN γ	IL-1 α	IL-1 β
		IL-2	IL-3	IL-4	IL-5
		IL-6	IL-10	IL-12p70	IL-17a
		MCP-1	MIP-1 α	RANTES	TNF α

Product Range

Analyte	Assay Type	Calibrator Range	Upper Limit of Quantification (ULOQ)	Lower Limit of Quantification (LLOQ)	Limit of Detection	Precision* (Inter-assay)	Precision* (Intra-assay)	Average* Linearity
GMCSF	Sandwich	560 - 0.77 (pg/mL)	392 (pg/mL)	1 (pg/mL)	0.38 (pg/mL)	4.9%	7.5%	94%
IFN γ	Sandwich	516 - 0.71 (pg/mL)	361 (pg/mL)	0.94 (pg/mL)	0.35 (pg/mL)	3.3%	5.5%	110%
IL-1 α	Sandwich	610 - 0.84 (pg/mL)	427 (pg/mL)	1.1 (pg/mL)	0.42 (pg/mL)	3.3%	7.4%	108%
IL-1 β	Sandwich	2500 - 3.43 (pg/mL)	1750 (pg/mL)	4.5 (pg/mL)	1.71 (pg/mL)	3.1%	9.3%	81%
IL-2	Sandwich	360 - 0.49 (pg/mL)	252 (pg/mL)	0.6 (pg/mL)	0.25 (pg/mL)	5.1%	6.6%	88%
IL-3	Sandwich	60 - 0.08 (pg/mL)	42 (pg/mL)	0.1 (pg/mL)	0.04 (pg/mL)	4.7%	7.0%	106%
IL-4	Sandwich	400 - 0.55 (pg/mL)	280 (pg/mL)	0.7 (pg/mL)	0.27 (pg/mL)	6.2%	6.6%	98%
IL-5	Sandwich	546 - 0.75 (pg/mL)	382 (pg/mL)	0.96 (pg/mL)	0.37 (pg/mL)	6.6%	5.2%	103%
IL-6	Sandwich	504 - 0.69 (pg/mL)	353 (pg/mL)	0.91 (pg/mL)	0.35 (pg/mL)	7.3%	6.1%	106%
IL-10	Sandwich	667 - 0.91 (pg/mL)	467 (pg/mL)	1.33 (pg/mL)	0.46 (pg/mL)	7.5%	3.0%	108%
IL-12p70	Sandwich	572 - 0.78 (pg/mL)	400 (pg/mL)	1.01 (pg/mL)	0.39 (pg/mL)	5.5%	10.8%	107%
IL-17a	Sandwich	760 - 1.04 (pg/mL)	532 (pg/mL)	1.4 (pg/mL)	0.52 (pg/mL)	6.5%	6.1%	92%
MCP-1	Sandwich	9400 - 12.89 (pg/mL)	6580 (pg/mL)	16.8 (pg/mL)	6.45 (pg/mL)	5.5%	9.2%	96%
MIP-1 α	Sandwich	1000 - 1.37 (pg/mL)	700 (pg/mL)	1.8 (pg/mL)	0.69 (pg/mL)	10.9%	6.7%	86%
RANTES	Sandwich	1653 - 2.27 (pg/mL)	1157 (pg/mL)	2.92 (pg/mL)	1.13 (pg/mL)	5.0%	5.5%	108%
TNF α	Sandwich	429 - 0.59 (pg/mL)	300 (pg/mL)	0.69 (pg/mL)	0.29 (pg/mL)	6.9%	5.2%	104%

* Actual values may vary from kit to kit. Please see the In-Kit Certificate of Analysis included in your kit for specific values.

The Q-Plex™ Array

Spotting in 96 Well Plate

- ▶ Robotic liquid handlers print 20-50nL spots of capture antibody
- ▶ Each spot is a unique assay within the well
- ▶ Low spot-to-spot variability (CV)
- ▶ Spot size 350-500µm
- ▶ Plates are QC'd for spot quality

Performing the Assay

- ▶ Add as little as 25µL of sample
- ▶ Wash
- ▶ Add mix of detection antibodies specific to kit
- ▶ Wash
- ▶ Add streptavidin conjugated HRP

Detection of Sample

- ▶ With the addition of substrate, a response is produced
- ▶ If antigen is present the spot emits a signal proportional to the amount of antigen in the sample
- ▶ If no antigen is present, the spot is not visible

Image Capture





- ▶ An image of the plate is taken via a high resolution camera contained within the Q-View™ Imager LS

Image Analysis

- ▶ Image is opened in Q-View™ Software
- ▶ Spots are automatically found on plate image
- ▶ Intensity of spot response is measured and raw data is generated
- ▶ User imports product specs and well layout

Data Analysis

- ▶ Raw data is analyzed and compared to in-plate standard
- ▶ Regression models used to calculate unknowns
- ▶ Standard curves are calculated and sample and statistical data is exported

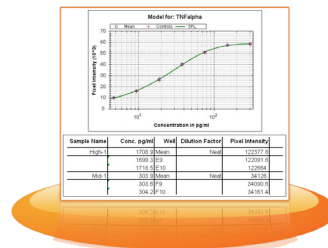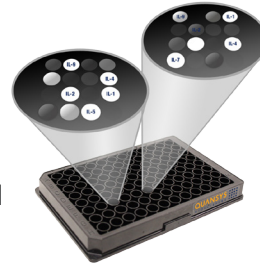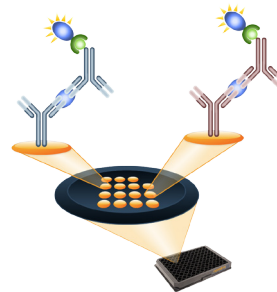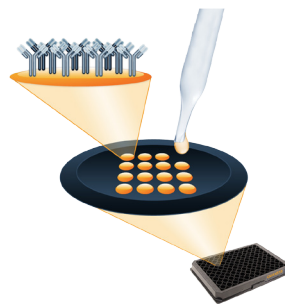

Increase productivity with Quansys products and services. Save time, sample, and money.

Q-Plex™ ARRAYS

Our Q-Plex Arrays are quantitative multiplex ELISAs with distinct proteins deposited in a defined array. Choose one of our standard kits for immediate delivery or customize the exact array you need.

Q-View™ IMAGERS

The Q-View™ Imager LS is a high quality, low-cost imaging systems for chemiluminescent assay imaging.

Q-View™ SOFTWARE

The Q-View Software is a user-friendly image analysis app that enables the acquisition and analysis of large amounts of multiplex ELISA data.

SAMPLE TESTING SERVICE

Our lab runs immunoassays on a wide variety of biological samples. Increase productivity and focus on your core research while we test your samples for you.



Quansys Biosciences
365 North 600 West
Logan, UT 84321
www.quansysbio.com

Fax: 435-750-6869
Phone: 435-752-0531
Toll Free: 888-QUANSYS (782-6797)

Quansys is an ISO 9001:2015 and ISO 13485:2016 registered company and complies with GMP. Products are designed, developed, and manufactured according to the procedures outlined in our Quality Management System.

©2022 Quansys Biosciences Inc. All rights reserved. Q-View and Q-Plex are Trademarks of Quansys Biosciences Inc. ISO is a trademark of the International Standards Organization. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.

QUANSYS
BIOSCIENCES



www.quansysbio.com