QUANSYS BIOSCIENCES

Q-Plex™ Human Innate Immunity (7-Plex)

Version 1.1
For Research Use Only
Not For Use in Diagnostic Procedures

TABLE OF CONTENTS

Name and Intended Use	2
Principle of the Assay	3
Sample Collection and Storage	4
Warnings and Precautions	5
Kit Contents & Storage	6
Materials Required but not Supplied	7
Abbreviated Protocol	8
Preparation	8
Running the Assay	8
Assay Preparation	9
Assay Procedure	10
Analyzing a Q-Plex™ Image	13
Appendix A: Plate Washing Methods	14
Automated Wash Method	14
Multichannel Pipette Wash Method	15
Plate Diagram	16

Symbol	Explanation
REF	Catalog Number
LOT	Lot Number
	Use By YYYY-MM-DD
	Temperature Limitation
	Manufacturer
	Keep Away from Sunlight

Manufactured and Distributed by:
Quansys Biosciences
365 North 600 West Logan, Utah 84321, USA
TEL: (888) 782-6797
E-MAIL: info@quansysbio.com

For International Distributors see:
http://www.quansysbio.com/distributors

QUANSYS BIOSCIENCES

NAME AND INTENDED USE

Q-Plex™Human Innate Immunity (7-Plex)
Quansys Biosciences Catalog Number 111233HU

The Q-Plex Human Innate Immunity (7-Plex) is a chemiluminescent Enzyme-Linked Immunosorbent Assay (ELISA) intended for quantitative concurrent measurement of GMCSF, IFNα,IFNß, IFNγ, IL-1α, IL-1ß, and IL-18 in serum and EDTA/Heparin Plasma samples.

PRINCIPLE OF THE ASSAY

This multiplex assay is based on the microplate sandwich enzyme immunoassay technique for the measurement of GMCSF, IFNα,IFNß, IFNγ, IL-1α, IL-1ß, and IL-18.

The GMCSF, IFNα,IFNß, IFNγ, IL-1α, IL-1ß, and IL-18 assays use 2 different antibodies specific for their respective targets. Samples or calibrators are pipetted into wells of a 96-well plate arrayed with immobilized, analyte-specific antibodies that capture GMCSF, IFNα,IFNß, IFNγ, IL-1α, IL-1ß, and IL-18. After washing away any unbound protein, a mixture that contains biotinylated analyte-specific antibodies is added. The biotinylated antibodies complete the sandwich for each specific arrayed analyte. After washing away unbound biotinylated antibody, streptavidin-horseradish peroxidase (SHRP) is added. Following an additional wash, the amount of SHRP bound to each location of the array is proportional to the amount of GMCSF, IFNα,IFNß, IFNγ, IL-1α, IL-1ß, and IL-18 initially captured.

The amount of conjugated enzyme on each location of the array is measured with the addition of a chemiluminescent substrate.

SAMPLE COLLECTION AND STORAGE

The sample collection and storage conditions are intended as general guidelines. The Clinical and Laboratory Standards Institute (CLSI GP44-A4) recommends the following storage conditions for samples: Samples should be stored at room temperature no longer than 8 hours. If the assay will not be completed within 8 hours, the samples should be refrigerated at +2°C to +8°C. If the assay will not be completed within 48 hours, or if the samples will be stored beyond 48 hours, samples should be frozen at -20°C or lower. Samples should not be repeatedly frozen and thawed. Frozen samples must be mixed well after thawing and prior to testing. Do not use bacterially contaminated samples. It is the responsibility of the individual laboratory to use all available references and/or its own studies to determine its own specific stability criteria.

Cell Culture Supernatant – Remove particulates by centrifugation and assay immediately or aliquot and store samples at ≤ -20°C. Avoid repeated freeze-thaw cycles.

Serum – Use a serum separator tube (SST) and allow samples to clot for 30 minutes before centrifugation for 15 minutes at 1000 x g. Remove serum and assay immediately or aliquot and store samples at ≤ -20°C. Avoid repeated freeze-thaw cycles.

Plasma – Collect plasma using EDTA, or Heparin as an anticoagulant. Centrifuge for 15 minutes at 1000 x g within 30 minutes of collection. Assay immediately or aliquot and store samples at ≤ -20°C. Avoid repeated freeze-thaw cycles.

WARNINGS AND PRECAUTIONS

1. Read all instructions before beginning test.
2. For research use only. Not for use in diagnostic procedures.
3. The kit should not be used beyond the expiration date on the kit label.
4. If running multiple kits, a calibration curve must be included for each 96-well plate. The curve from one plate cannot be used to calculate sample values from other plates.
5. Do not mix or substitute reagents with those from other kits or lots.
6. Pre-wet pipette tips 3 times by drawing up the liquid into the pipette and then dispensing back into the original vessel prior to the addition of samples, calibrators, or external controls to the microplate.
7. Load all calibrators, controls, and samples into the microplate within 10 minutes of each other.
8. Be exact with incubation times, particularly the streptavidin-horseradish peroxidase (SHRP) incubation.
9. Be exact when mixing Substrate A and B+ and mix thoroughly.
10. Do not allow the substrate or SHRP to be exposed to UV light, as this may degrade these reagents.
11. Do not allow the plate to dry out between steps.
12. All products are carefully validated; however due to the variability encountered in biological sample matrices, the possibility of interference or sample matrix effects cannot be excluded.

    Warning: The calibrator contains components of human origin. These components have been tested at the donor level and found to be negative for HBsAg, HIV-1 and HIV-2 antibodies, and HCV. (Or these components have been heat inactivated.) However, consider all materials as potentially infectious and use only approved guidelines for the proper handling and disposal of infectious material.

KIT CONTENTS & STORAGE

Unopened Kit - Store at 2-8°C. Do not use past kit expiration date. Do not re-use.

Part/Description	Storage of opened/reconstituted material
Q-Plex™ Array Microplate REF 111252HU 1 arrayed and blocked 96-well polystyrene micro-titer plate	2-8°C until kit expiration
Wash Buffer Concentrate (20X) REF 101158GR Liquid, 50 mL/vial of a concentrated solution of buffered surfactant	
Sample Diluent REF 101104HU Liquid, 10 mL/vial of a buffered protein solution with heterophilic antibody and rheumatoid factor blockers and preservatives	
Detection Mix REF (111256HU) Liquid, 6 mL/vial of biotinylated antibodies in a buffered protein solution with preservatives	
Calibrator 1 REF 102655HU Lyophilized, recombinant antigens in a buffered protein base	2-8°C until kit expiration Discard unused reconstituted calibrator.
Calibrator 2 REF 102855HU Lyophilized, recombinant antigens in a buffered protein base	2-8°C until kit expiration Discard unused reconstituted calibrator.
Streptavidin-HRP REF 101173GR Liquid, 6 mL/vial of streptavidin-conjugated horse-radish peroxidase	Do not expose to UV light. 2-8°C until kit expiration
Substrate A REF 101159GR Liquid, 3.5 mL/vial	Do not expose to UV light. Store mixed substrate solution at room temperature (20-25°C) for up to 1 week. Store unmixed solutions at 2-8°C until kit expiration.
Substrate B+ REF 101120GR Liquid, 3.5 mL/vial	
Plate Seals (3) REF S1064 Adhesive strips	Non-perishable

MATERIALS REQUIRED BUT NOT SUPPLIED

In addition to the kit contents listed, the following materials are required to run this assay:

1. Multichannel pipette (20-200 µL) and/or single channel pipettes (20-1000 µL) with appropriate tips.
2. 10 mL serological pipette.
3. Polypropylene tubes or polypropylene 96-well plate(s) for sample and calibrator preparation.
    a. Example: Nunc® MicroWell™ 96-Well Plates, Polypropylene, 249944; Eppendorf® LoBind Protein or Genomic Microcentrifuge Tubes, 022431102.
4. Q-View™ Imager and Q-View™ Software.
5. Microplate shaker capable of 500 RPM.
6. Deionized water.
7. Microplate washer (An optional multichannel pipette protocol is provided in Appendix A).
8. 1 L graduated cylinder for the preparation of wash buffer.

ABBREVIATED PROTOCOL

Preparation

1. Install Q-View Software (page 9).
2. Set up the imager (page 9).
3. Set up microplate washer (page 14) and shaker (page 9).
4. Reconstitute and prepare reagents (page 9).

Running the Assay

1. Prepare the calibration curve using the prepared Calibrator and Sample Diluent according to the In-Kit Certificate of Analysis (page 10).
2. Prepare the samples with Sample Diluent (page 10).
3. Load the calibration curve and samples onto the plate. Shake for 120 minutes (500 RPM) at room temperature (page 11).
4. Wash the plate 3 times, add the Detection Mix, and shake for 60 minutes (500 RPM) at room temperature (page 11).
5. Wash the plate 3 times, add the Streptavidin-HRP, and shake for 20 minutes (500 RPM) at room temperature (page 11).
6. Wash the plate 6 times and add the prepared Substrate (page 11).
7. Capture and analyze image of the plate (page 11).

ASSAY PREPARATION

1. Install Q-View Software on any computers that will be used for analysis or operating a Q-View Imager.
2. Set up the imager. For imager-specific instructions, see www.quansysbio.com/manuals.
3. Set the microplate shaker to 500 RPM.
4. Allow all kit components to come to room temperature (20-25°C) prior to use.
5. Prepare Wash Buffer: Place 50 mL of Wash Buffer Concentrate (20X) into 950 mL deionized water and mix thoroughly.
6. Prepare Calibrator: Reconstitute lyophilized Calibrator 1 using Sample Diluent with the volume specified on the In-Kit Certificate of accompanying the kit. Allow Calibrator 1 to sit for 5 minutes and mix thoroughly. Reconstitute Calibrator 2 using Calibrator 1 with the volume specified on the In-Kit Certificate of Analysis. Allow Calibrator 2 to sit for 5 minutes. Mix thoroughly. Use prepared Calibrator within 30 minutes.
7. Prepare chemiluminescent substrate: Combine 3 mL of Substrate A with 3 mL of Substrate B+ and mix gently. **Do not expose to UV light. Store at room temperature (20-25°C) after mixing.**

ASSAY PROCEDURE


Allow all reagents to equilibrate to room temperature (20-25°C) before use and prepare as directed in the previous sections. It is recommended that all calibrators, samples, and controls be assayed in duplicate.

1. Using the previously prepared Calibrator (see "Assay Preparation" step 6), prepare an 8-point calibration curve (7 points plus 1 blank) in either polypropylene tubes or a polypropylene 96-well plate (see "Materials Required but Not Supplied" for example tubes and plates).

    a. Pipette 200 µL of prepared Calibrator into the first tube or well.

    b. Pipette 120 µL of Sample Diluent into the other 7 tubes or wells.

    c. Transfer 60 µL of undiluted prepared Calibrator from the first tube or well into the second, mix thoroughly, and repeat the transfer from tube to tube or well to well for 5 more points, leaving the last tube or well without any prepared Calibrator. This process is diagrammed below. The undiluted prepared Calibrator serves as the high point of the calibration curve. The Sample Diluent serves as the blank.

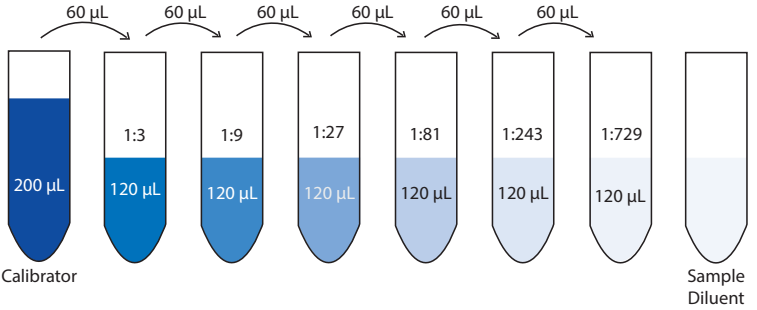

2. Prepare samples by diluting 1:2. 1 part sample to 1 part Sample Diluent) in either polypropylene tubes or a polypropylene 96-well plate in a volume sufficient to provide 50 µL of diluted sample per well.

Note: If you anticipate that your analyte concentrations will be greater than the high point of the calibration curves as reported on the In-Kit Certificate of Analysis, use the sample diluent to dilute your samples further.

3. Add 50 µL per well of the prepared calibration curve, controls, and samples to the Q-Plex™ Array Microplate. Load all samples and calibration curve to the plate within 10 minutes.
4. Cover the plate with a plate seal and incubate on a plate shaker set to 500 RPM for 120 minutes at room temperature (20-25˚C).
5. Wash the plate 3 times (see Appendix A).
6. Add 50 µL per well of Detection Mix, cover with a new plate seal, and return to the plate shaker set to 500 RPM for 60 minutes at room temperature (20-25˚C).
7. Wash the plate 3 times (see Appendix A).
8. Add 50 µL per well of Streptavidin-HRP, cover with a new plate seal, and return to the plate shaker set to 500 RPM for 20 minutes at room temperature (20-25°C).
9. Wash the plate 6 times (see Appendix A).
10. Add 50 µL per well of previously prepared substrate (see "Assay Preparation" step 7). Image plate immediately. Wait no longer than 5 minutes to commence imaging.
11. Place the plate in the Q-View Imager.
12. Open Q-View Software, create or open a project, and click Acquire Image.
13. Use the chart below to determine the recommended imaging time based on imager model:

Imager Model	Recommended Imaging Time
Q-View™ Imager LS	270 seconds
Q-View™ Imager Pro/Plus	300 seconds

14. Click the Capture Image(s) button. Users may perform Well Assignment while images are being captured.

    Note: Details about these imaging steps are available in the Q-View Software Manual viewable at www.quansysbio.com/manuals or within Q-View Software under **Support > Manual**.

15. Dispose of all used and unused materials in accordance with local regulations.

ANALYZING A Q-PLEX™ IMAGE


The following summarizes a general workflow for analyzing a Q-Plex image in Q-View Software. Each of these steps is described in greater detail in the Q-View Software and Imager manuals, viewable at www.quansysbio.com/support or within Q-View Software under **Support > Manual**.

1. Acquire or import an image into Q-View Software as previously described.
2. Enter the **Software Product Code** (found on the In-Kit Certificate of Analysis) into the **Product Code** field of the software.
3. **Image Processing:** Align the plate overlay as follows:
    a. Set the overlay: If using the **Auto-Set Plate Overlay** feature, this will occur automatically. Otherwise, go to **Overlay Options > Set Plate Overlay**.
    b. To visualize bright or dim spots, optimize the display using **Image Options > Adjust Gamma** (does not affect the data).
    c. Optimize overlay alignment: Go to **Overlay Options > Adjust Plate Overlay** to pivot the overlay, **Adjust Well** and **Adjust Spot** to move individual wells and spots, then **Auto-Adjust Spots** to automatically snap each circle of the overlay to the nearest spot of the image beneath.
4. **Well Assignment:** Label wells as samples, calibrators, controls, or negatives (blank), and specify their dilution factors. Use Templates to quickly assign layouts that are repeated often or export the layout as a .csv file.
5. **Data Analysis:** Once **Image Processing** and **Well Assignment** are complete, select **Data Analysis**. Click **Perform Analysis** to generate charts, concentrations, and statistics based on optimal regression settings provided by the product definition.

APPENDIX A: PLATE WASHING METHODS

Automated Wash Method

1. Use a program that will aspirate and dispense 300-400 µL wash buffer.
2. Ensure that the plate washer has a program that will not scratch the bottom of the microplate but will prevent plate drying. This is critical to prevent damage to the capture antibody arrays. The simplest method to avoid plate drying is to leave a small, uniform amount (1-3 µL) of wash in the well after the final aspiration and add the next reagent to the plate as quickly as possible.

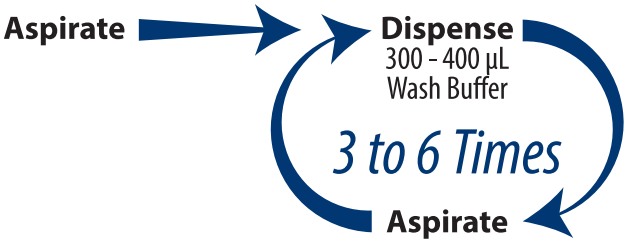

3. Leaving a uniform amount of wash in the wells can be accomplished with an automated washer by slightly increasing the aspiration height of the washer head.

For example:

Process	Distance	Steps on a Biotek ELX-405
Aspiration Height	3.81 mm	30
Aspiration Position	1.28 mm from center	-28
Dispense Height	15.24 mm	120

No soak or shake cycles are needed during wash steps


4. Connect the prepared wash buffer to your automatic plate washer.
5. Run 1-2 priming cycles to make sure that the wash buffer is running through the plate washer. When the buffer has run through the machine, the waste will be foamy.
6. In a spare microtiter plate, dispense 100 µL wash buffer, ensure that all pins dispense uniformly, then aspirate and ensure that all pins aspirate uniformly. This will ensure that all pins are functioning.
7. Prime the plate washer 1 time before the first wash step. When running the assay, perform the wash 3 or 6 times according to the protocol.

Multichannel Pipette Wash Method

1. Just prior to washing, pour the prepared wash buffer into a trough or tray.
2. After each incubation, flick the solution out of the plate over a waste container before starting the wash protocol.
3. Using a multichannel pipette, dispense 400 µL of wash buffer into each of the wells used in the test.
4. Aggressively flick the wash buffer out over a waste container.
5. This washes the plate 1 time. When the assay procedure calls for 3 or 6 washes, repeat steps 3-4 accordingly.
6. Tap the plate upside down on a paper towel to remove any residual wash.
7. Proceed immediately to dispense the next solution so drying does not occur.

We take great care to ensure that customers have success using our products and services. If you have further questions about running the assay, data analysis, troubleshooting, or our products or services, please contact us at 1-888-QUANSYS (1-888-782-6797) or at support@quansysbio.com.

PLATE DIAGRAM

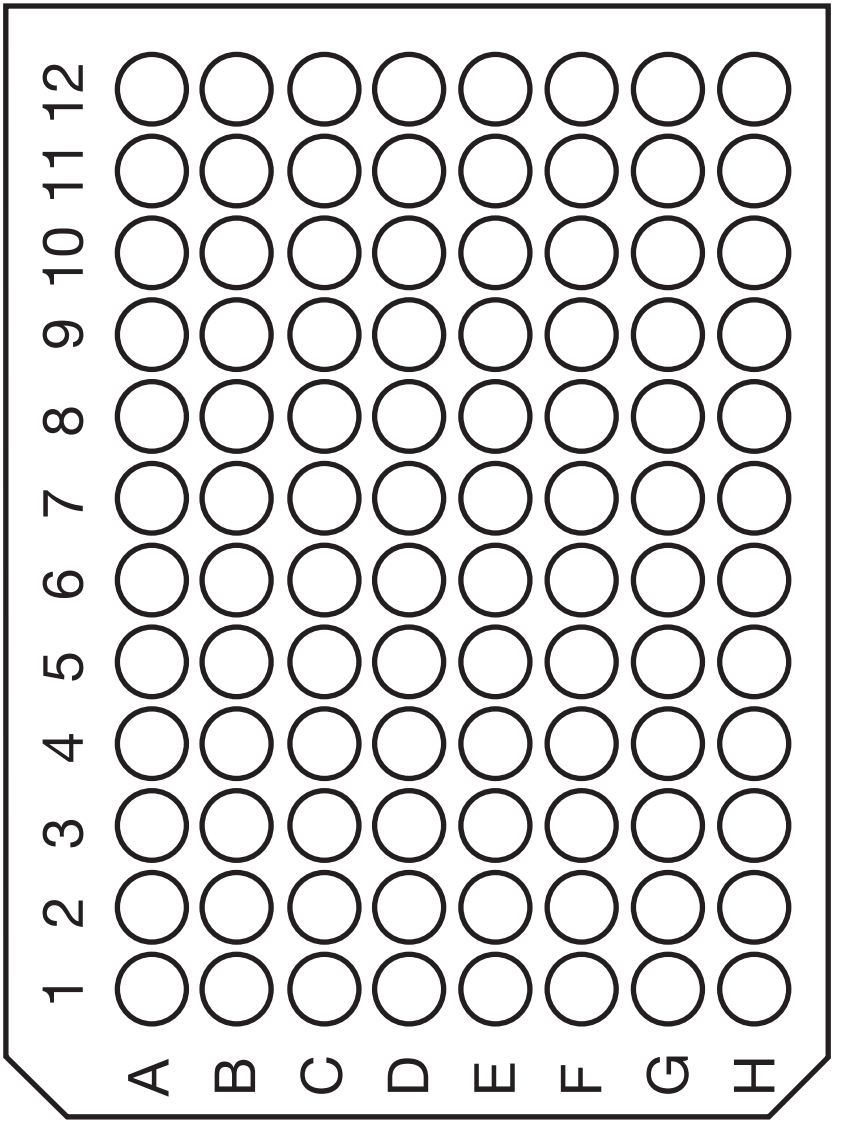

For Research Use Only

Q-Plex™ ARRAY

For technical support, email: support@quansysbio.com

Part Number 111217HU

365 North 600 West, Logan, Utah 84321
Tel: 1-888-782-6797 | Fax: (435)750-6869 | www.quansysbio.com